

Digital Microscope Simple Operation Guide

Digital microscope simple operation guide

1. Open the package and take out the microscope, and connect the microscope to the computer host (make sure the microscope and the computer socket are plugged in place without looseness)

2. Install the software: Please log in to <http://soft.hvscam.com> to download the software and select the HiViewSetup.exe icon to open it, then click Next---I Accept---"Install---"Finish

3. Double-click the desktop icon to open the software. After entering the software, you will see the following interface:

1. Open the package and take out the microscope, and connect the microscope to the computer mainframe (make sure that the sockets of the microscope and computer are in place and not loose)



2. Install the software: Please login to <http://soft.hvscam.com> to download the software, select and open the icon HiViewSetup.exe , and then click Next---Receive --- Install ---Finish

3. Double click the icon HiViewSetup.exe on the desktop to open the software and after entering the software, you will see the following interface:



4. Click "Device" in the menu bar and then click "Open" and select the microscope device name, color, compression, and image resolution, etc.:


4. Click "Device" in the menu bar and then click "Open" and select the microscope device name, color/compression, and image resolution, and so on:

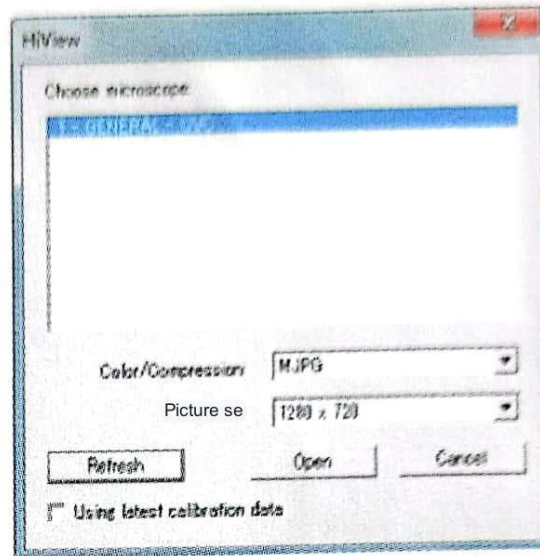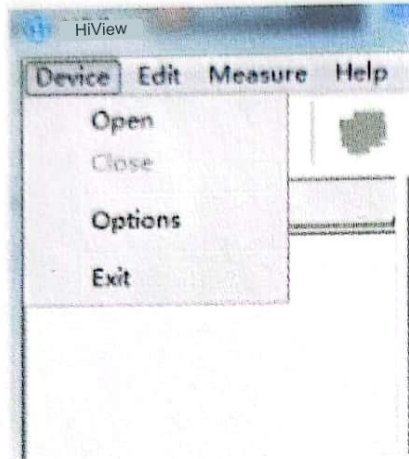

5. Click "Open" to enter the system page to perform various operations, as shown in the following figure: (Caution! Microscope and software details


For instructions, please download the electronic manual).

5. Click "Open" to enter the system page to implement various operations, as shown below: (Note! Please download the manual for detailed instructions on the microscope and software)

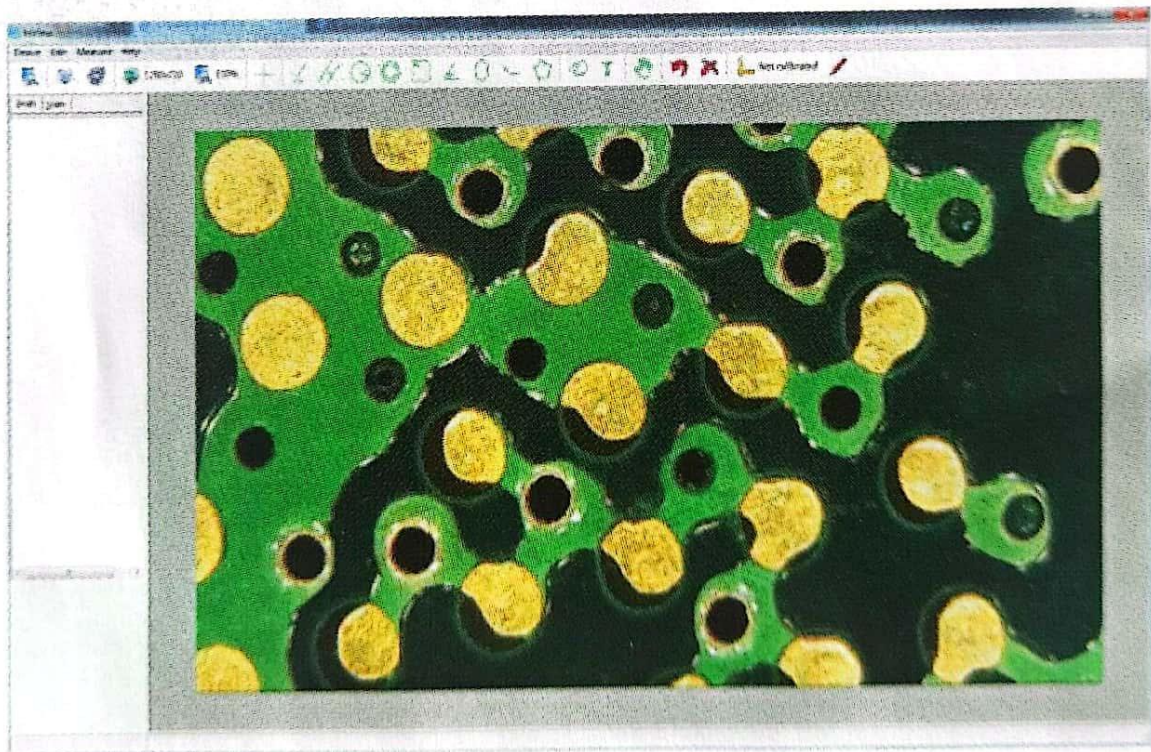

6. Measuring and line

drawing and other operational examples are shown below:

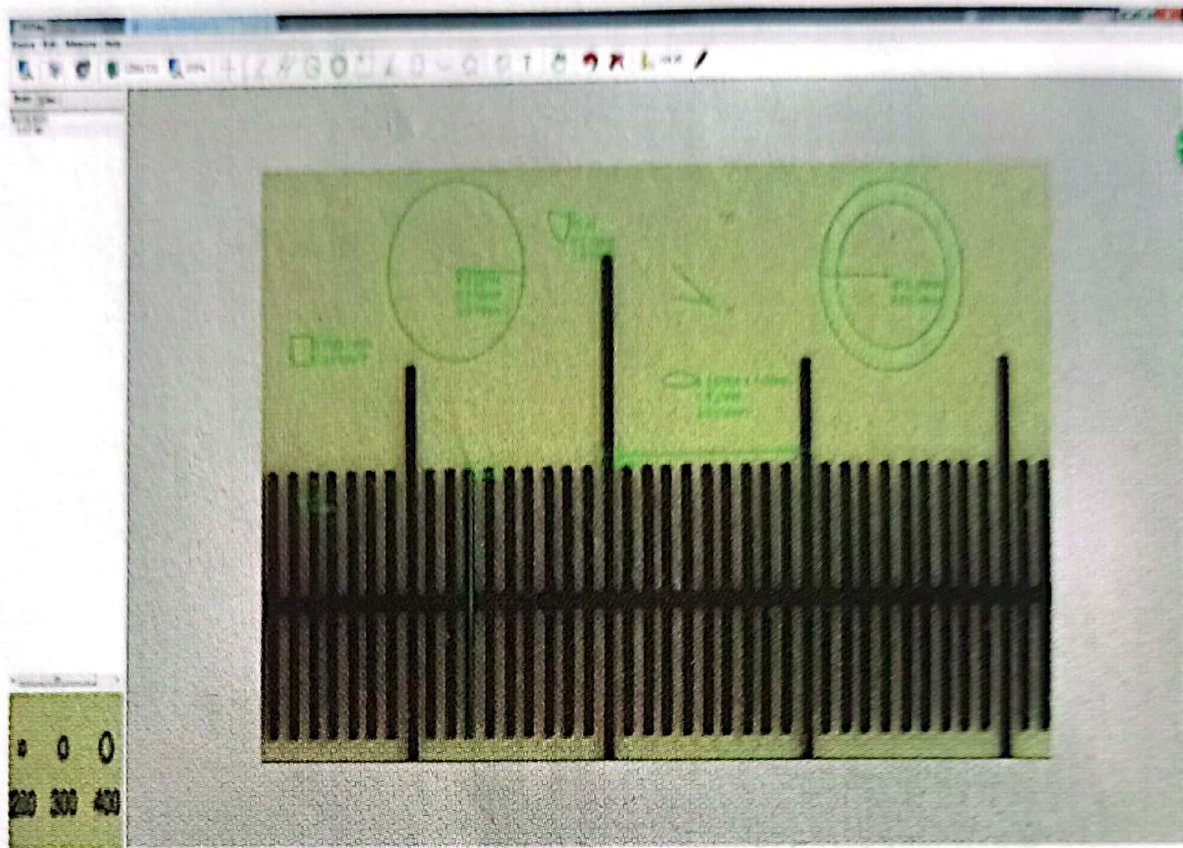

Note: Due to product software upgrades, the ZOOM button on the original body cannot be used. Please use this function at the "Modify Display Scale" in the software. Note: Due to product software upgrade, the "ZOOM" button on the body cannot be used. Please use "ZOOM" in the software.